



## Run Info

Host Name	GXB01275 (localhost)
Position	X2
Experiment Name	kokia
Sample ID	Kc1
Run ID	12343096-8afe-4a8a-9833-0d6e1fd161a1
Acquisition ID(s)	a138d8742ac26d27c894fa4e50a41b62e53e8976, 5f7b1a6d8732f8dc18aef21562c09b1134bfdbc6
Flow Cell Id	FAT12726
Start Time	April 1, 15:07
Run Length	8d 20h 56m

## Run Summary

Reads Generated	2.4 M
Passed Bases	10.07 Gb
Failed Bases	2.08 Gb
Estimated Bases	13.47 Gb
Percentage Basecalled	102%


## Run Parameters

Flow Cell Type	FLO-MIN112
Kit	SQK-LSK112
Initial bias voltage	-200 mV
FAST5 output	Enabled
FASTQ output	Enabled
BAM output	Disabled
Bulk file output	Disabled
Active channel selection	Enabled
Basecalling	Enabled
Specified run length	72 hours
FAST5 reads per file	4000
FAST5 output options	vbz_compress,fastq,raw
FASTQ reads per file	4000
FASTQ output options	compress
Mux scan period	1 hour 30 minutes
Reserved pores	0 %
Basecall model	dna_r10.4_e8.1_sup.cfg
Read filtering	min_qscore=10
Read splitting	enable=on


## Versions

MinKNOW	21.11.7
MinKNOW Core	4.5.4
Bream	6.3.5
Guppy	5.1.13

Cumulative Output Reads

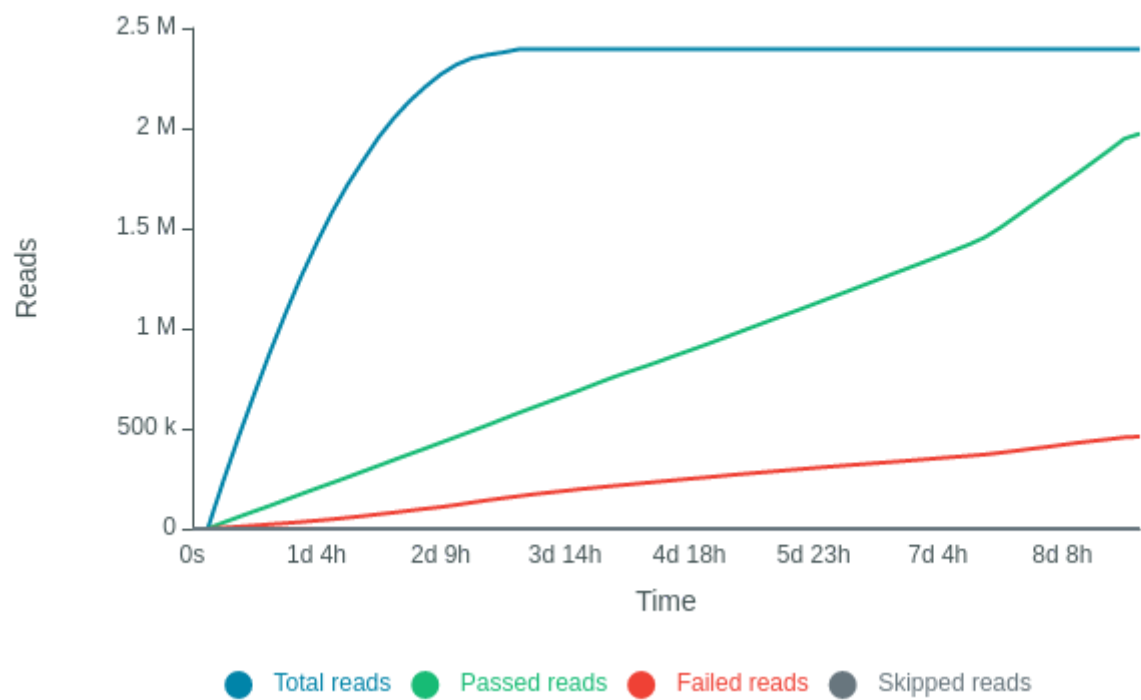


Cumulative Output Bases

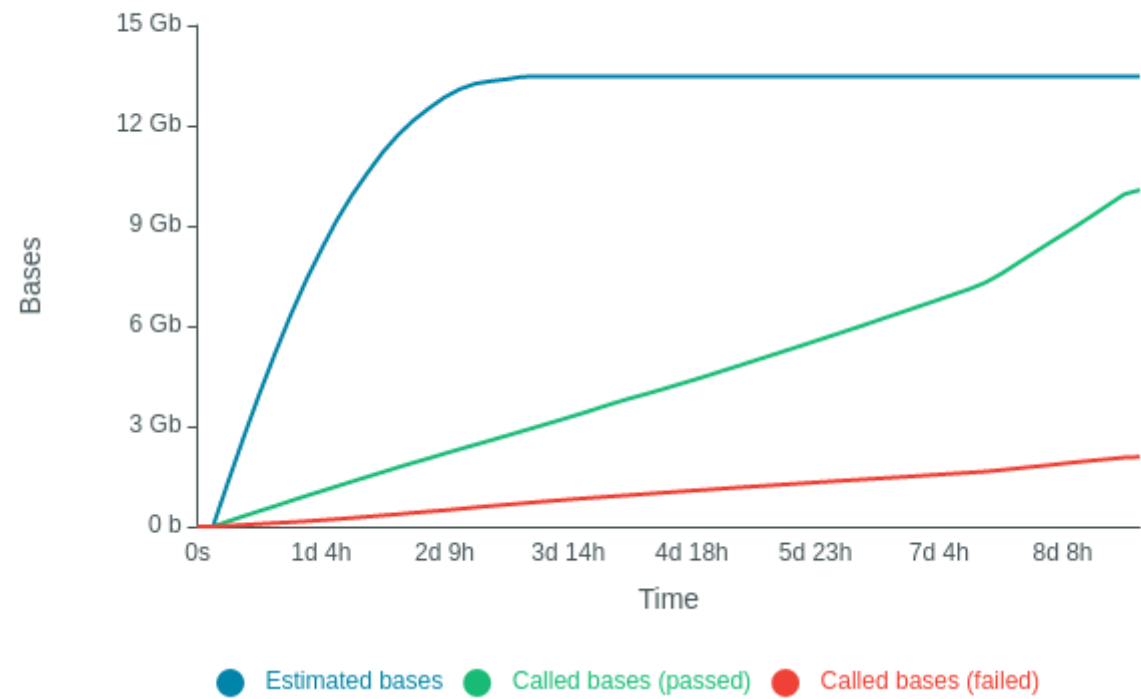


**Read Length Histogram Estimated Bases - Outliers Discarded**

Estimated N50: 8.35 kb

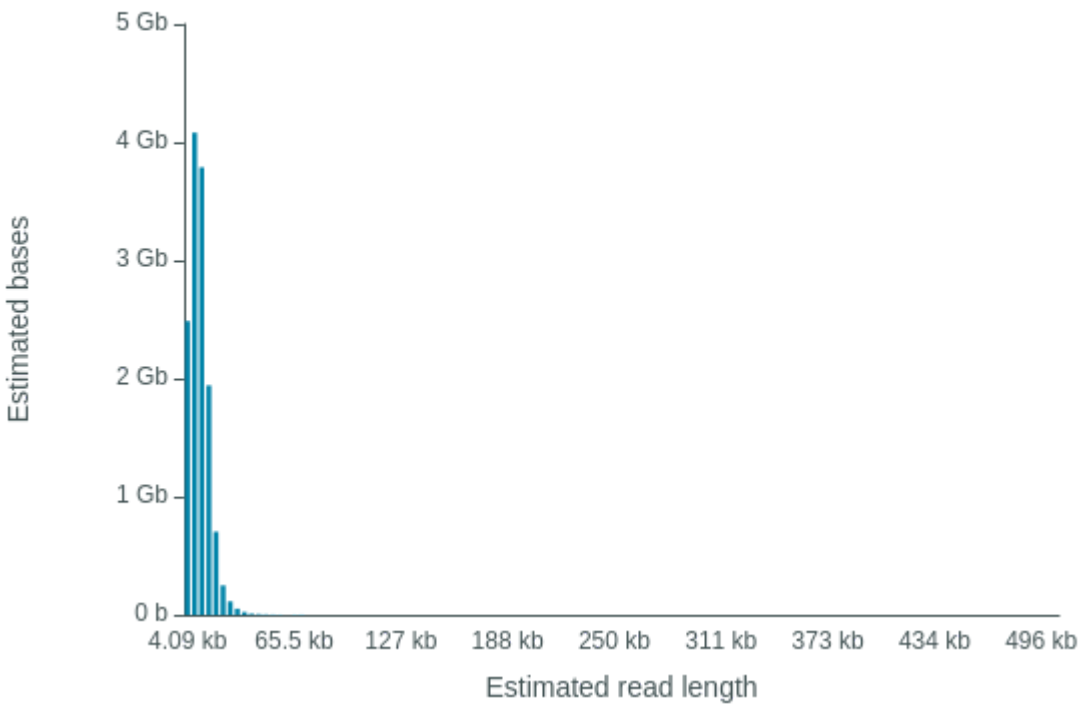


**Read Length Histogram Basecalled Bases - Outliers Discarded**

Estimated N50: 7.51 kb

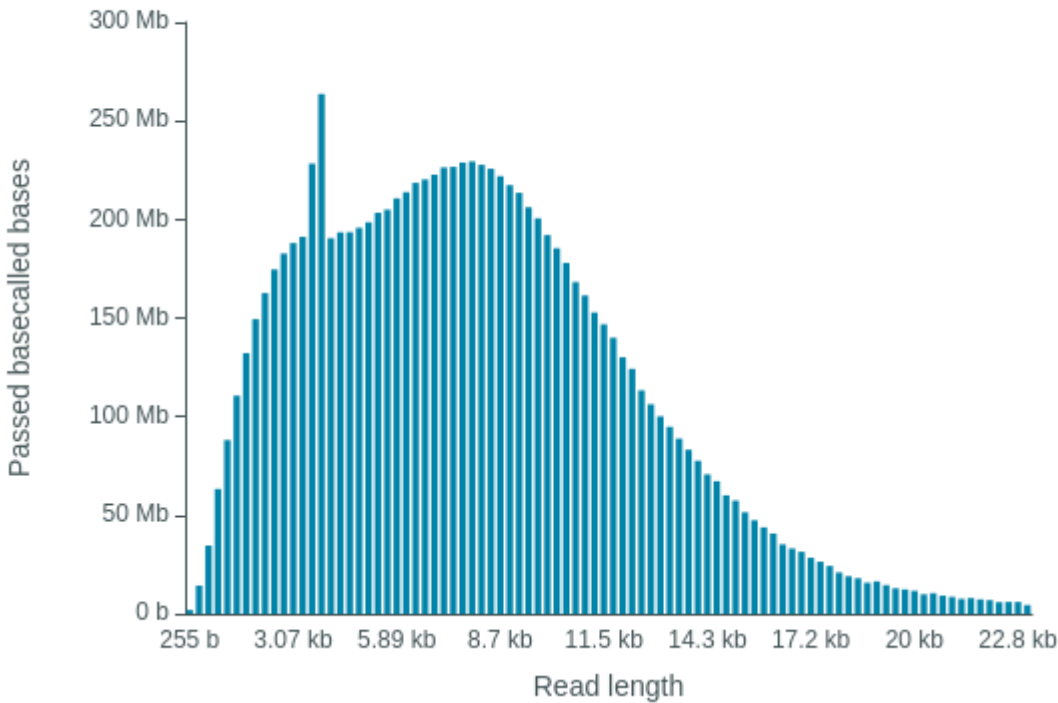


**Read Length Histogram Estimated Bases**

Estimated N50: 8.35 kb

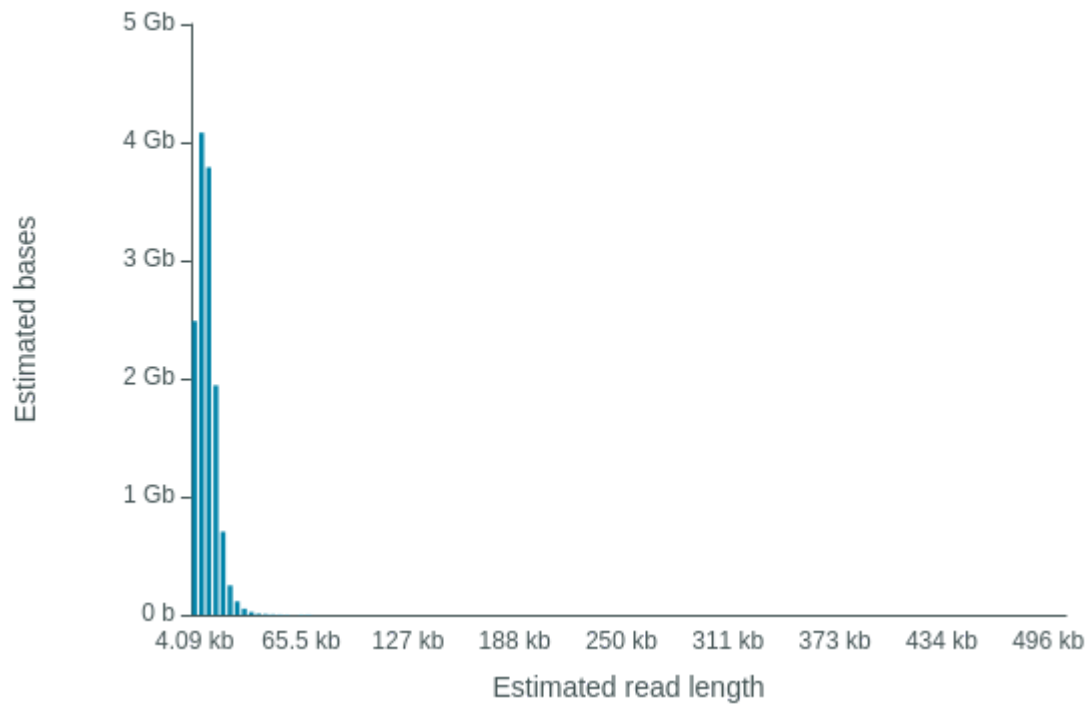


**Read Length Histogram Basecalled Bases**



Estimated N50: 7.51 kb

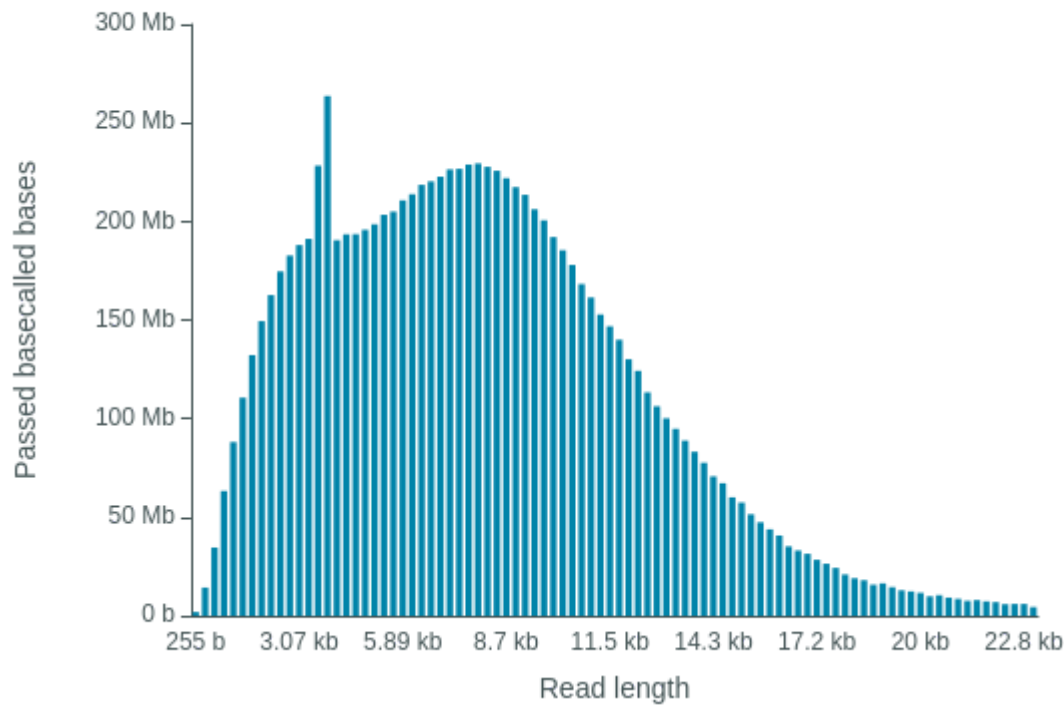


### Duty Time Grouped

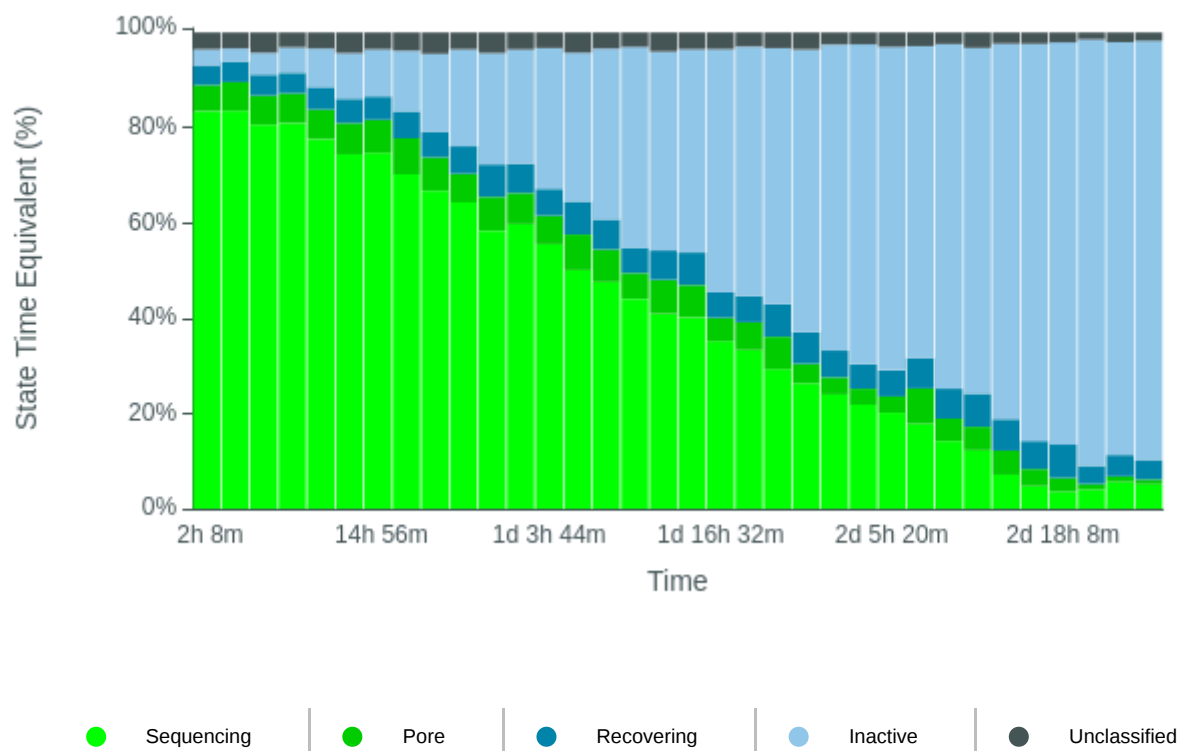


### Duty time Categorised

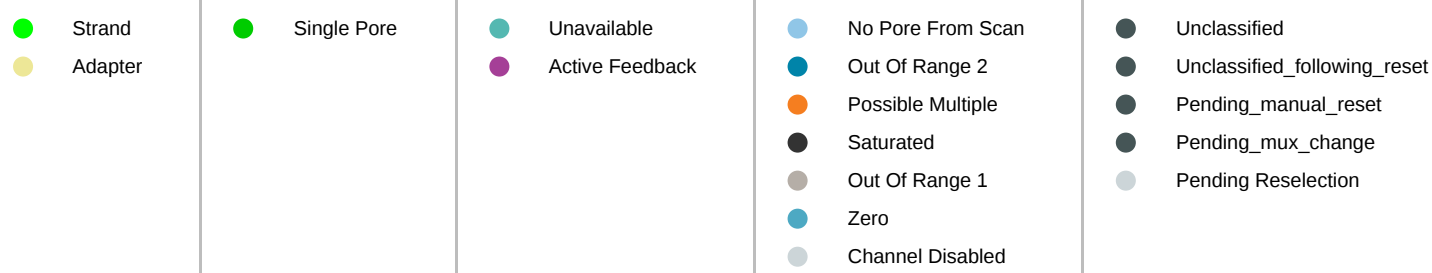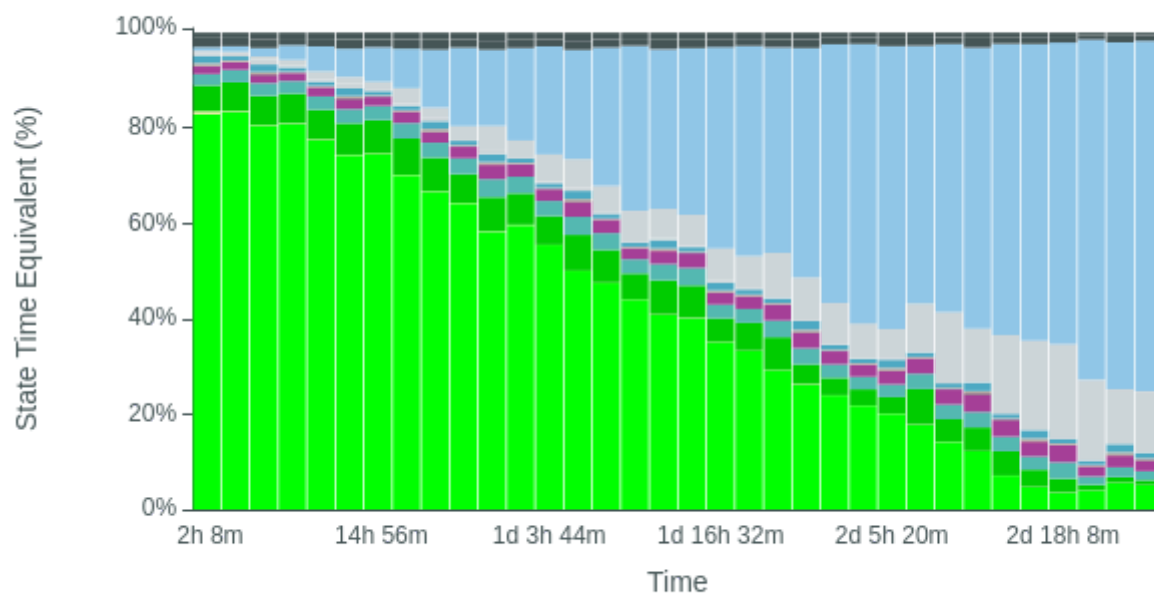


Mux Scan Grouped

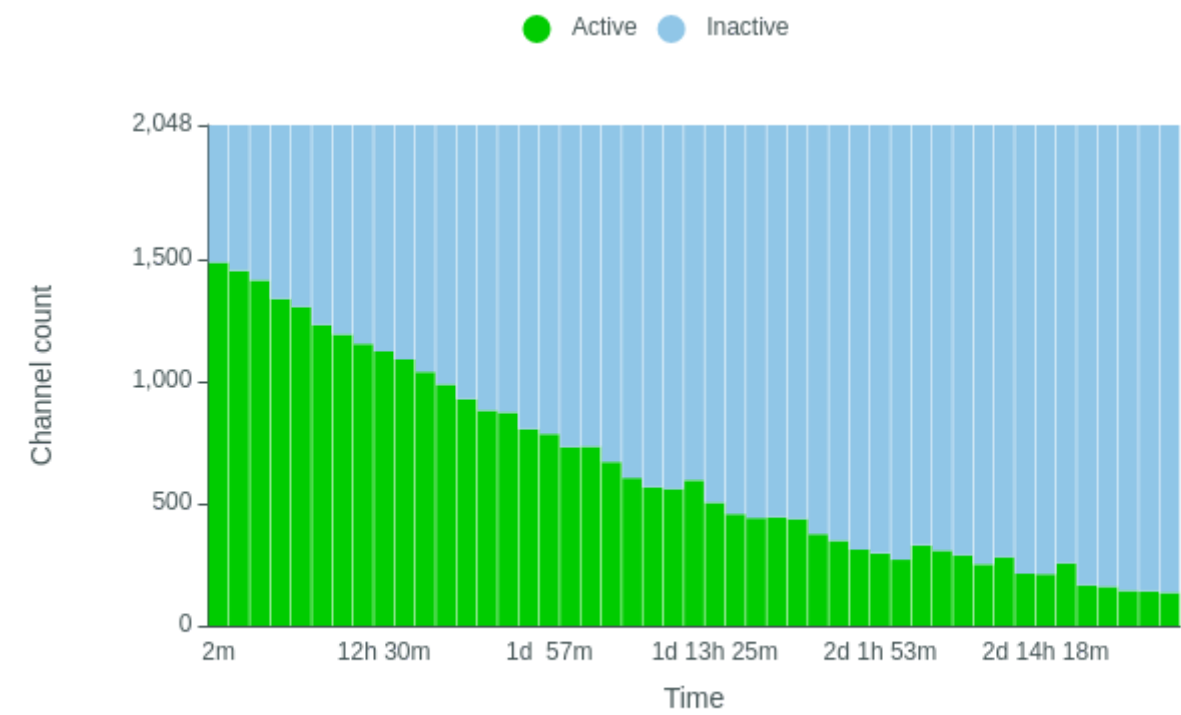


Mux Scan Categorised

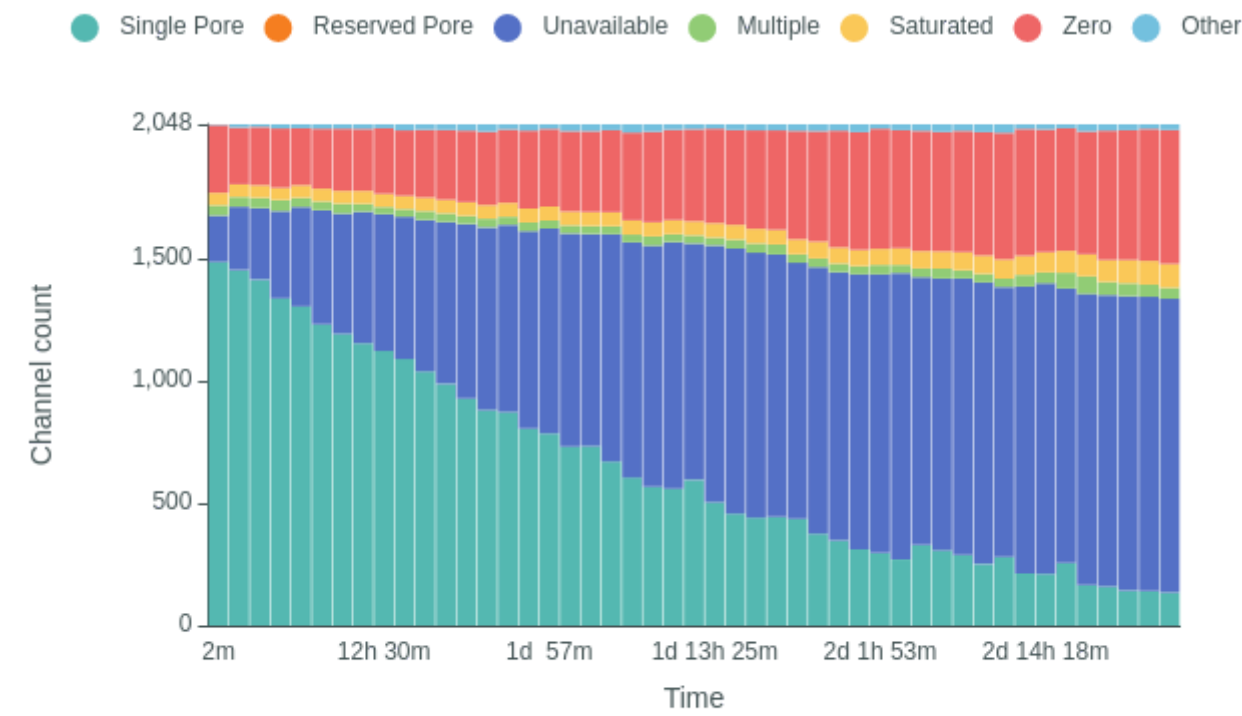


Temperature History.

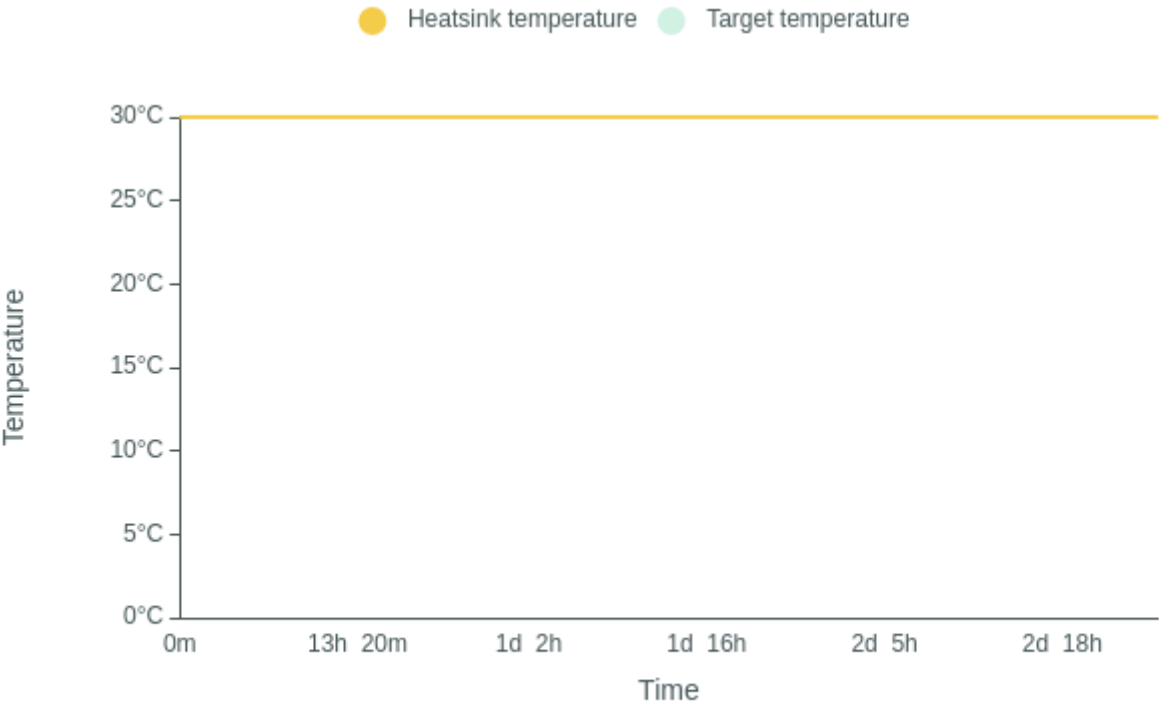


Bias Voltage History.

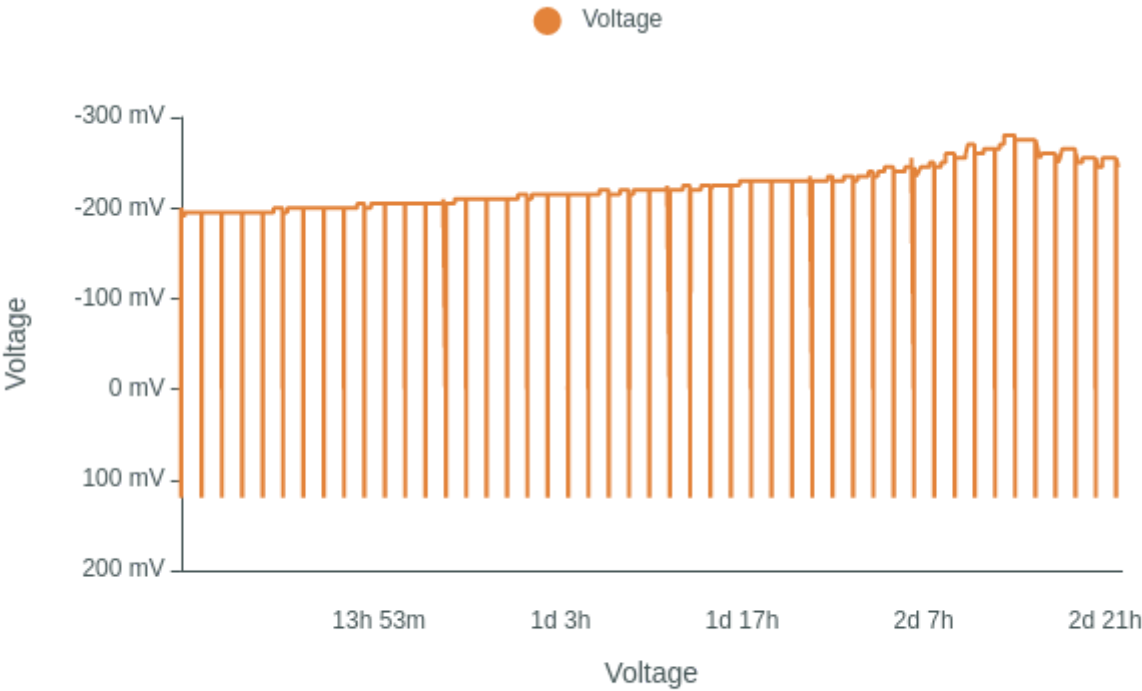

Translocation Speed

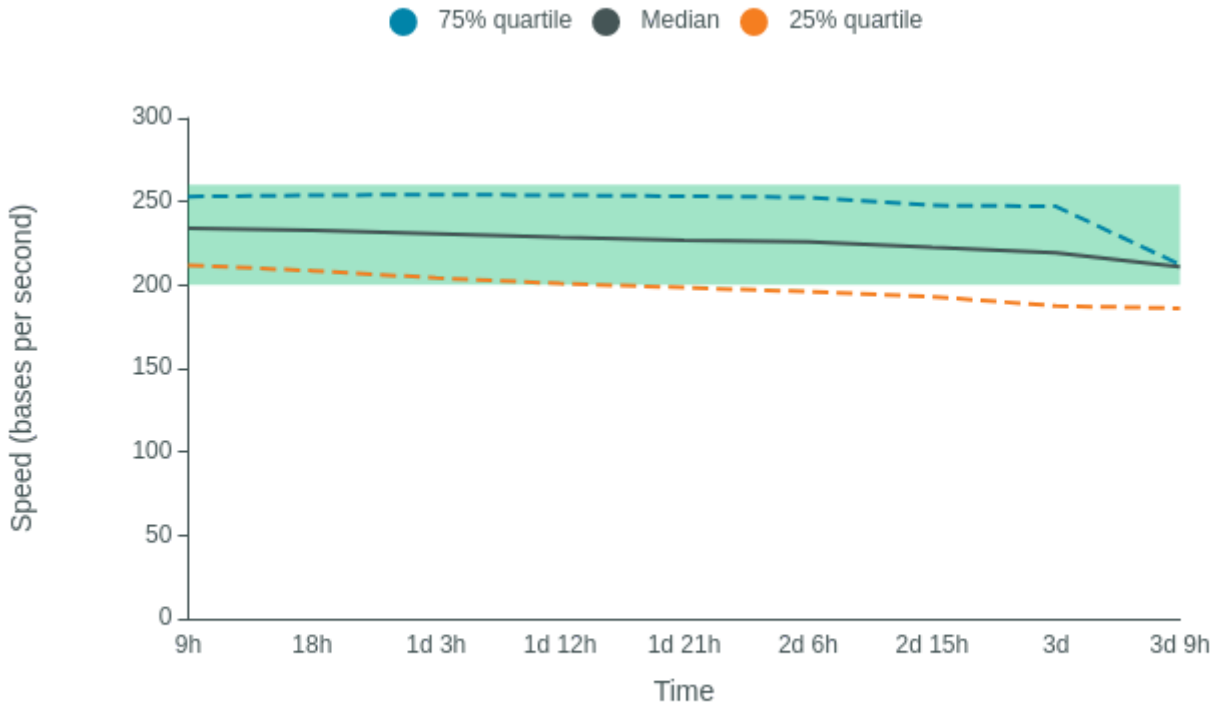


QScore

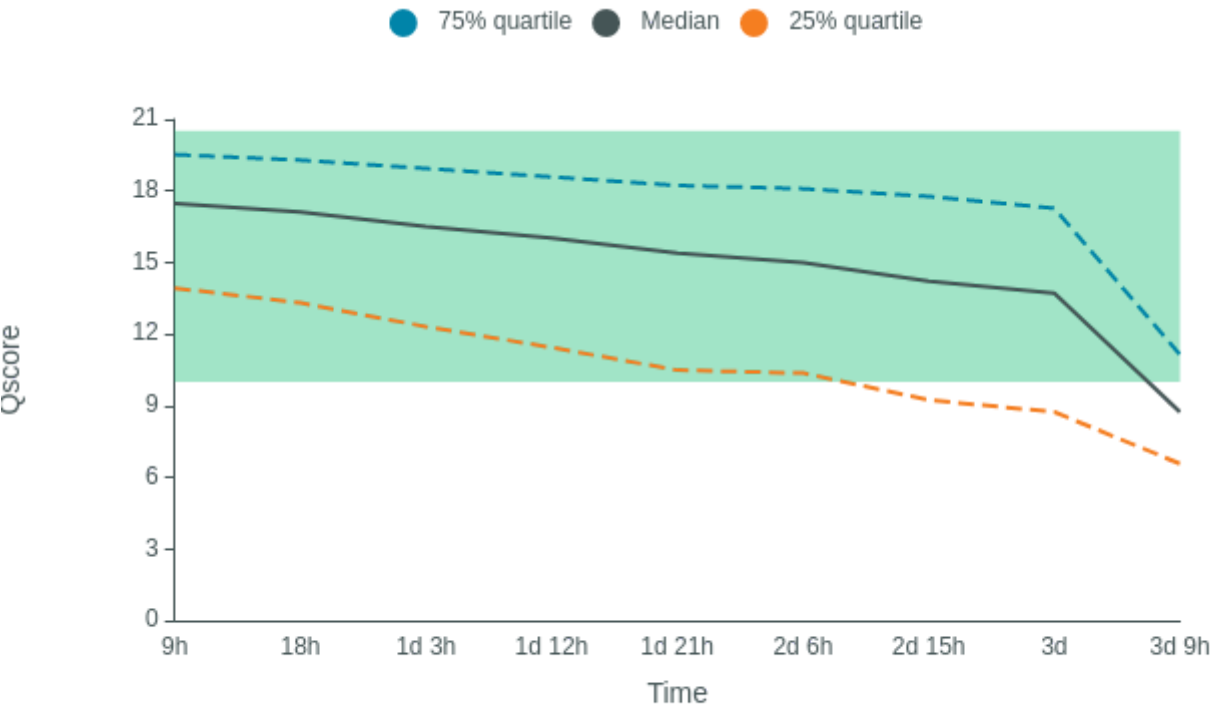



Disk Write Performance

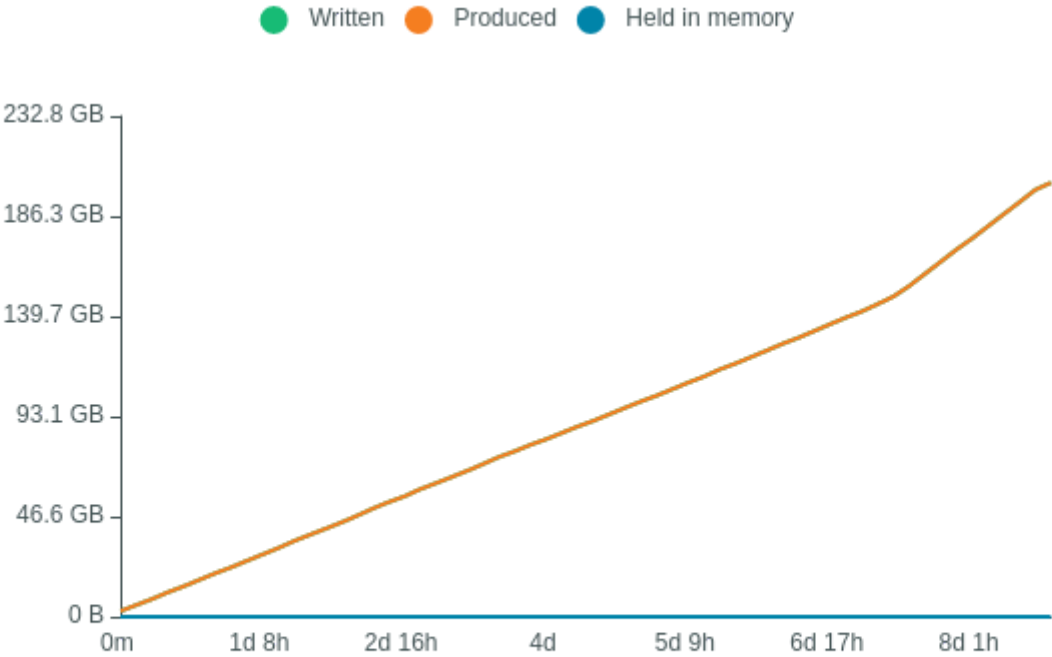

**Run Debug Messages**

- The sequencing run has finished, but basecalling may continue April 4, 15:12
- Mux scan for flow cell FAT12726 has found a total of 135 pores. 123 pores available for immediate sequencing April 4, 14:45
- Performing Mux Scan April 4, 14:42
- Mux scan for flow cell FAT12726 has found a total of 141 pores. 127 pores available for immediate sequencing April 4, 13:12
- Performing Mux Scan April 4, 13:10
- Mux scan for flow cell FAT12726 has found a total of 143 pores. 123 pores available for immediate sequencing April 4, 11:39
- Performing Mux Scan April 4, 11:37
- Mux scan for flow cell FAT12726 has found a total of 160 pores. 137 pores available for immediate sequencing April 4, 10:07
- Performing Mux Scan April 4, 10:05
- Mux scan for flow cell FAT12726 has found a total of 166 pores. 144 pores available for immediate sequencing April 4, 08:35
- Performing Mux Scan April 4, 08:32
- Mux scan for flow cell FAT12726 has found a total of 256 pores. 203 pores available for immediate sequencing April 4, 07:02
- Performing Mux Scan April 4, 07:00
- Mux scan for flow cell FAT12726 has found a total of 211 pores. 181 pores available for immediate sequencing April 4, 05:30
- Performing Mux Scan April 4, 05:27
- Mux scan for flow cell FAT12726 has found a total of 214 pores. 172 pores available for immediate sequencing April 4, 03:57
- Performing Mux Scan April 4, 03:55
- Mux scan for flow cell FAT12726 has found a total of 281 pores. 217 pores available for immediate sequencing April 4, 02:24
- Performing Mux Scan April 4, 02:22
- Mux scan for flow cell FAT12726 has found a total of 251 pores. 192 pores available for immediate sequencing April 4, 00:51
- Performing Mux Scan April 4, 00:48
- Mux scan for flow cell FAT12726 has found a total of 289 pores. 217 pores available for immediate sequencing April 3, 23:18
- Performing Mux Scan April 3, 23:15
- Mux scan for flow cell FAT12726 has found a total of 308 pores. 223 pores available for immediate sequencing April 3, 21:45
- Performing Mux Scan April 3, 21:42
- Mux scan for flow cell FAT12726 has found a total of 329 pores. 235 pores available for immediate sequencing April 3, 20:11
- Performing Mux Scan April 3, 20:09
- Mux scan for flow cell FAT12726 has found a total of 270 pores. 187 pores available for immediate sequencing April 3, 18:38
- Performing Mux Scan April 3, 18:35
- Mux scan for flow cell FAT12726 has found a total of 296 pores. 204 pores available for immediate sequencing April 3, 17:04
- Performing Mux Scan April 3, 17:02
- Mux scan for flow cell FAT12726 has found a total of 312 pores. 210 pores available for immediate sequencing April 3, 15:31
- Performing Mux Scan April 3, 15:28
- Mux scan for flow cell FAT12726 has found a total of 348 pores. 239 pores available for

- immediate sequencing April 3, 13:58
- Performing Mux Scan April 3, 13:55
- Mux scan for flow cell FAT12726 has found a total of 375 pores. 250 pores available for immediate sequencing April 3, 12:24
- Performing Mux Scan April 3, 12:22
- Mux scan for flow cell FAT12726 has found a total of 437 pores. 287 pores available for immediate sequencing April 3, 10:51
- Performing Mux Scan April 3, 10:48
- Mux scan for flow cell FAT12726 has found a total of 444 pores. 287 pores available for immediate sequencing April 3, 09:17
- Performing Mux Scan April 3, 09:15
- Mux scan for flow cell FAT12726 has found a total of 441 pores. 279 pores available for immediate sequencing April 3, 07:44
- Performing Mux Scan April 3, 07:41
- Mux scan for flow cell FAT12726 has found a total of 457 pores. 282 pores available for immediate sequencing April 3, 06:10
- Performing Mux Scan April 3, 06:08
- Mux scan for flow cell FAT12726 has found a total of 504 pores. 307 pores available for immediate sequencing April 3, 04:37
- Performing Mux Scan April 3, 04:34
- Mux scan for flow cell FAT12726 has found a total of 595 pores. 349 pores available for immediate sequencing April 3, 03:03
- Performing Mux Scan April 3, 03:01
- Mux scan for flow cell FAT12726 has found a total of 560 pores. 331 pores available for immediate sequencing April 3, 01:30
- Performing Mux Scan April 3, 01:27
- Mux scan for flow cell FAT12726 has found a total of 567 pores. 322 pores available for immediate sequencing April 2, 23:56
- Performing Mux Scan April 2, 23:54
- Mux scan for flow cell FAT12726 has found a total of 605 pores. 349 pores available for immediate sequencing April 2, 22:23
- Performing Mux Scan April 2, 22:20
- Mux scan for flow cell FAT12726 has found a total of 670 pores. 368 pores available for immediate sequencing April 2, 20:49
- Performing Mux Scan April 2, 20:47
- Mux scan for flow cell FAT12726 has found a total of 734 pores. 394 pores available for immediate sequencing April 2, 19:16
- Performing Mux Scan April 2, 19:13
- Mux scan for flow cell FAT12726 has found a total of 732 pores. 387 pores available for immediate sequencing April 2, 17:42
- Performing Mux Scan April 2, 17:40
- Mux scan for flow cell FAT12726 has found a total of 784 pores. 403 pores available for immediate sequencing April 2, 16:09
- Performing Mux Scan April 2, 16:06
- Mux scan for flow cell FAT12726 has found a total of 805 pores. 414 pores available for immediate sequencing April 2, 14:35
- Performing Mux Scan April 2, 14:33
- Mux scan for flow cell FAT12726 has found a total of 873 pores. 431 pores available for immediate sequencing April 2, 13:02
- Performing Mux Scan April 2, 13:00
- Mux scan for flow cell FAT12726 has found a total of 881 pores. 416 pores available for immediate sequencing April 2, 11:29

- Performing Mux Scan April 2, 11:26
- Mux scan for flow cell FAT12726 has found a total of 927 pores. 436 pores available for immediate sequencing April 2, 09:55
- Performing Mux Scan April 2, 09:53
- Mux scan for flow cell FAT12726 has found a total of 988 pores. 450 pores available for immediate sequencing April 2, 08:22
- Performing Mux Scan April 2, 08:19
- Mux scan for flow cell FAT12726 has found a total of 1038 pores. 467 pores available for immediate sequencing April 2, 06:48
- Performing Mux Scan April 2, 06:46
- Mux scan for flow cell FAT12726 has found a total of 1089 pores. 473 pores available for immediate sequencing April 2, 05:15
- Performing Mux Scan April 2, 05:12
- Mux scan for flow cell FAT12726 has found a total of 1122 pores. 475 pores available for immediate sequencing April 2, 03:41
- Performing Mux Scan April 2, 03:39
- Mux scan for flow cell FAT12726 has found a total of 1153 pores. 481 pores available for immediate sequencing April 2, 02:08
- Performing Mux Scan April 2, 02:05
- Mux scan for flow cell FAT12726 has found a total of 1192 pores. 484 pores available for immediate sequencing April 2, 00:34
- Performing Mux Scan April 2, 00:32
- Mux scan for flow cell FAT12726 has found a total of 1231 pores. 486 pores available for immediate sequencing April 1, 23:01
- Performing Mux Scan April 1, 22:58
- Mux scan for flow cell FAT12726 has found a total of 1304 pores. 501 pores available for immediate sequencing April 1, 21:27
- Performing Mux Scan April 1, 21:25
- Mux scan for flow cell FAT12726 has found a total of 1337 pores. 502 pores available for immediate sequencing April 1, 19:54
- Performing Mux Scan April 1, 19:51
- Mux scan for flow cell FAT12726 has found a total of 1414 pores. 505 pores available for immediate sequencing April 1, 18:20
- Performing Mux Scan April 1, 18:18
- Mux scan for flow cell FAT12726 has found a total of 1453 pores. 507 pores available for immediate sequencing April 1, 16:47
- Performing Mux Scan April 1, 16:44
- Mux scan for flow cell FAT12726 has found a total of 1486 pores. 508 pores available for immediate sequencing April 1, 15:13
- Performing Mux Scan April 1, 15:11
- Starting sequencing procedure April 1, 15:11
- Waiting up to 300 seconds for temperature to stabilise at 30.0°C April 1, 15:07