



## Run Info

Host Name	GXB01275 (localhost)
Position	X5
Experiment Name	kokia
Sample ID	Kk
Run ID	f1023326-45e9-4888-82fb-759d41a81fae
Acquisition ID(s)	d5c44a2f2ae97dc502e3a4fbe736f9c53f350f7f, 6636659f3a934d6f532d7b80b0c4959afca4a4f2
Flow Cell Id	FAU08590
Start Time	June 10, 14:15
Run Length	3d 0h 3m

## Run Summary

Reads Generated	1.88 M
Passed Bases	6.66 Gb
Failed Bases	1.85 Gb
Estimated Bases	9.86 Gb
Percentage Basecalled	101%


## Run Parameters

Flow Cell Type	FLO-MIN112
Kit	SQK-LSK112
Initial bias voltage	-200 mV
FAST5 output	Enabled
FASTQ output	Enabled
BAM output	Disabled
Bulk file output	Disabled
Active channel selection	Enabled
Basecalling	Enabled
Specified run length	72 hours
FAST5 reads per file	4000
FAST5 output options	vbz_compress,fastq,raw
FASTQ reads per file	4000
FASTQ output options	compress
Mux scan period	1 hour 30 minutes
Reserved pores	0 %
Basecall model	dna_r10.4_e8.1_hac.cfg
Read filtering	min_qscore=9
Read splitting	enable=on


## Versions

MinKNOW	21.11.7
MinKNOW Core	4.5.4
Bream	6.3.5
Guppy	5.1.13

Cumulative Output Reads

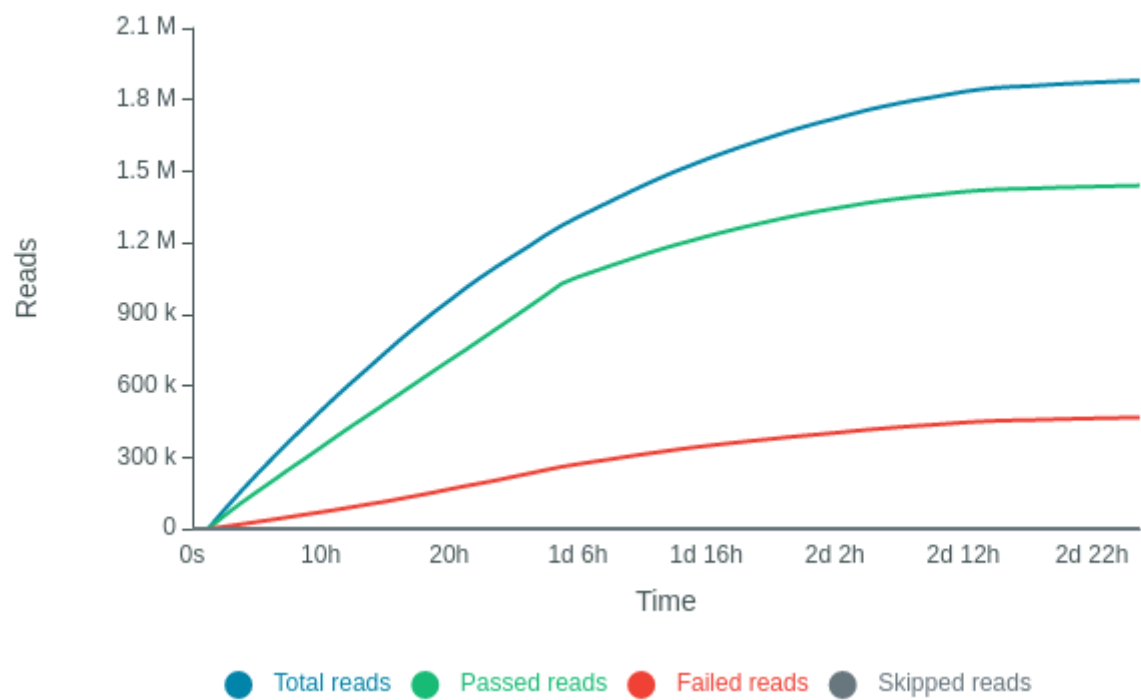


Cumulative Output Bases

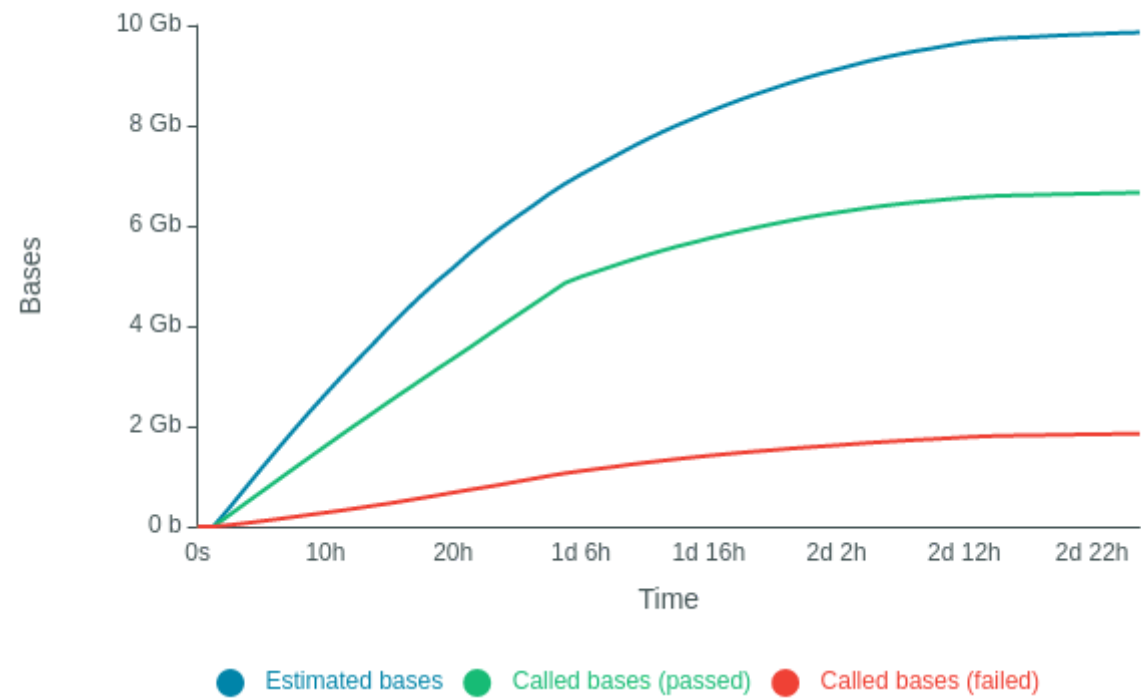


**Read Length Histogram Estimated Bases - Outliers Discarded**

Estimated N50: 8.3 kb

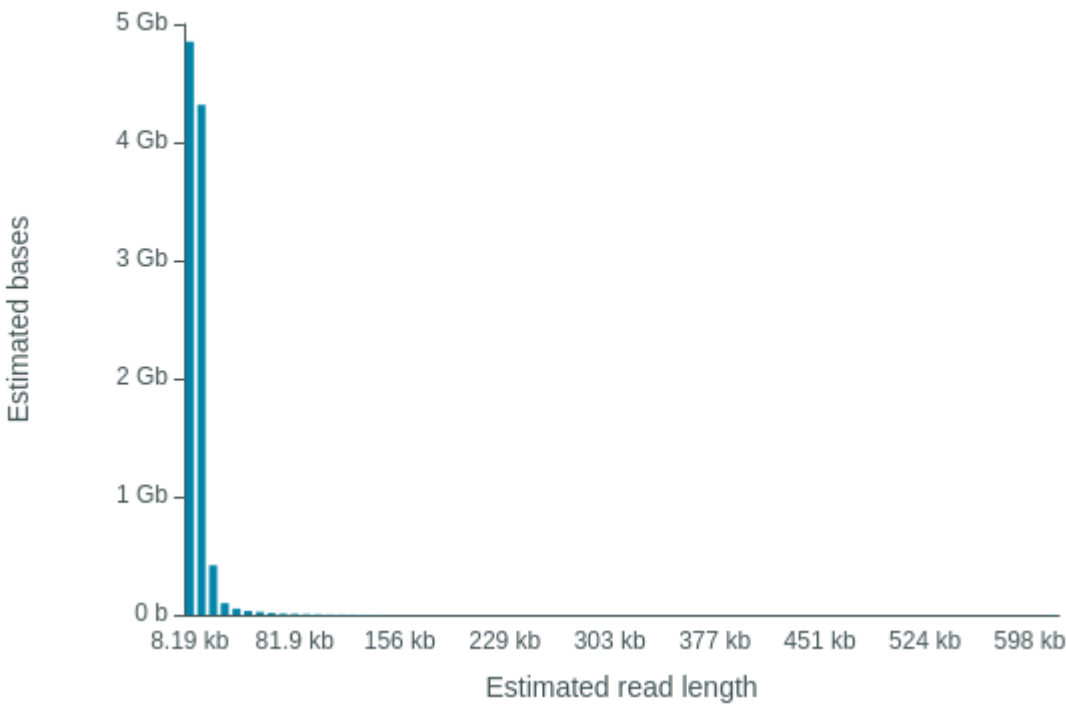


**Read Length Histogram Basecalled Bases - Outliers Discarded**

Estimated N50: 7.25 kb

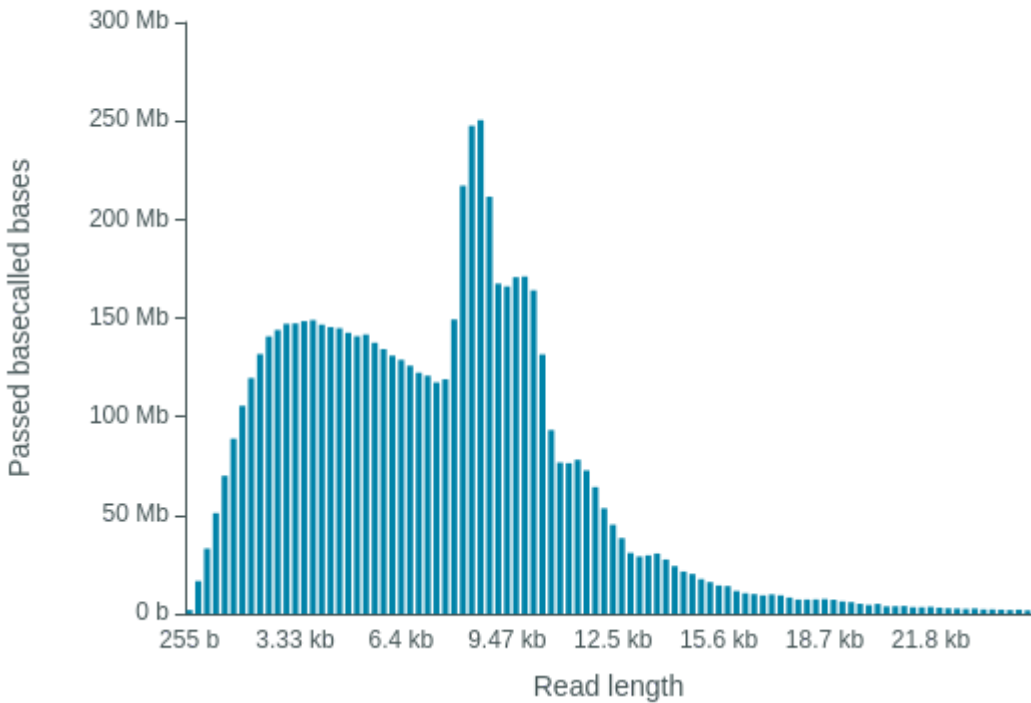


**Read Length Histogram Estimated Bases**

Estimated N50: 8.3 kb

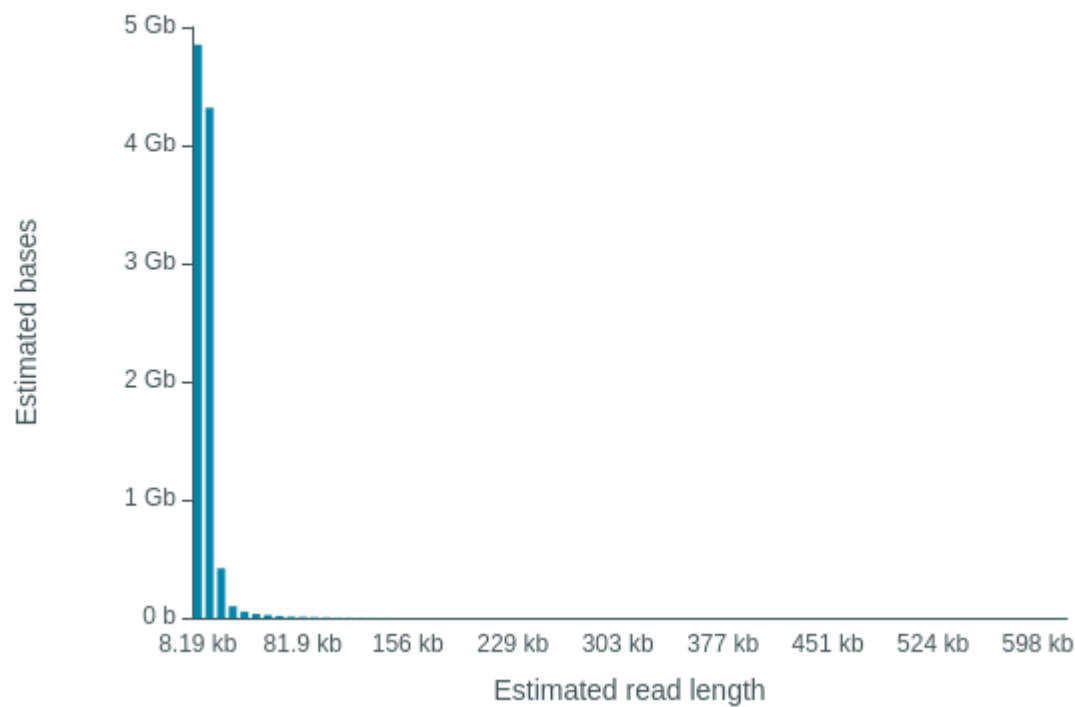


**Read Length Histogram Basecalled Bases**


Estimated N50: 7.25 kb

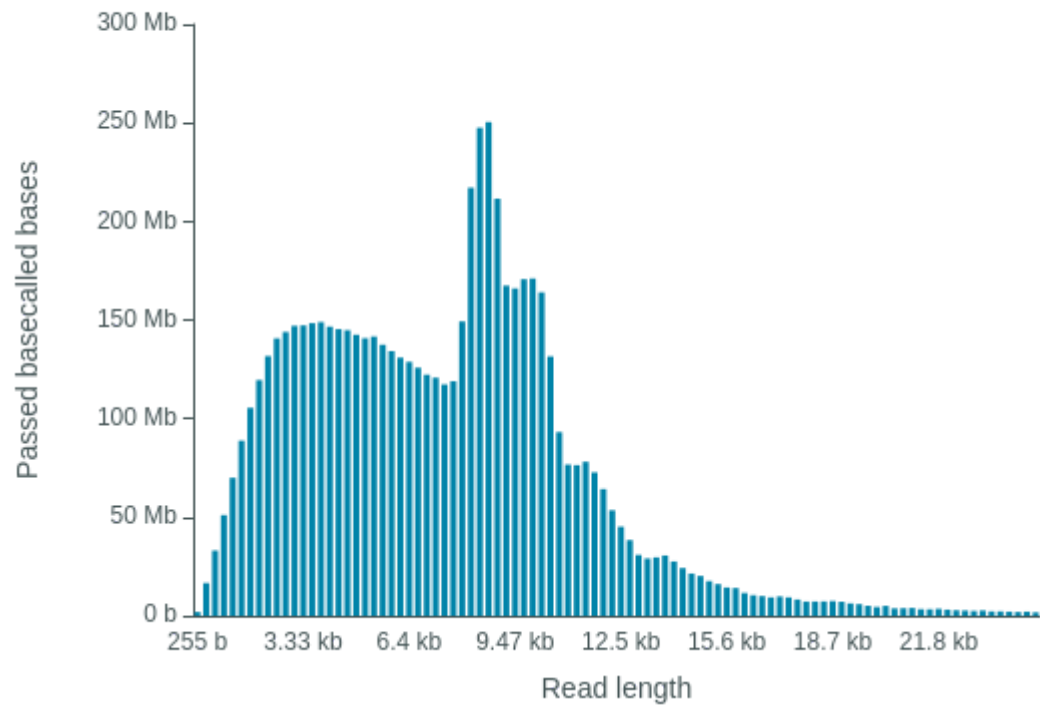


### Duty Time Grouped

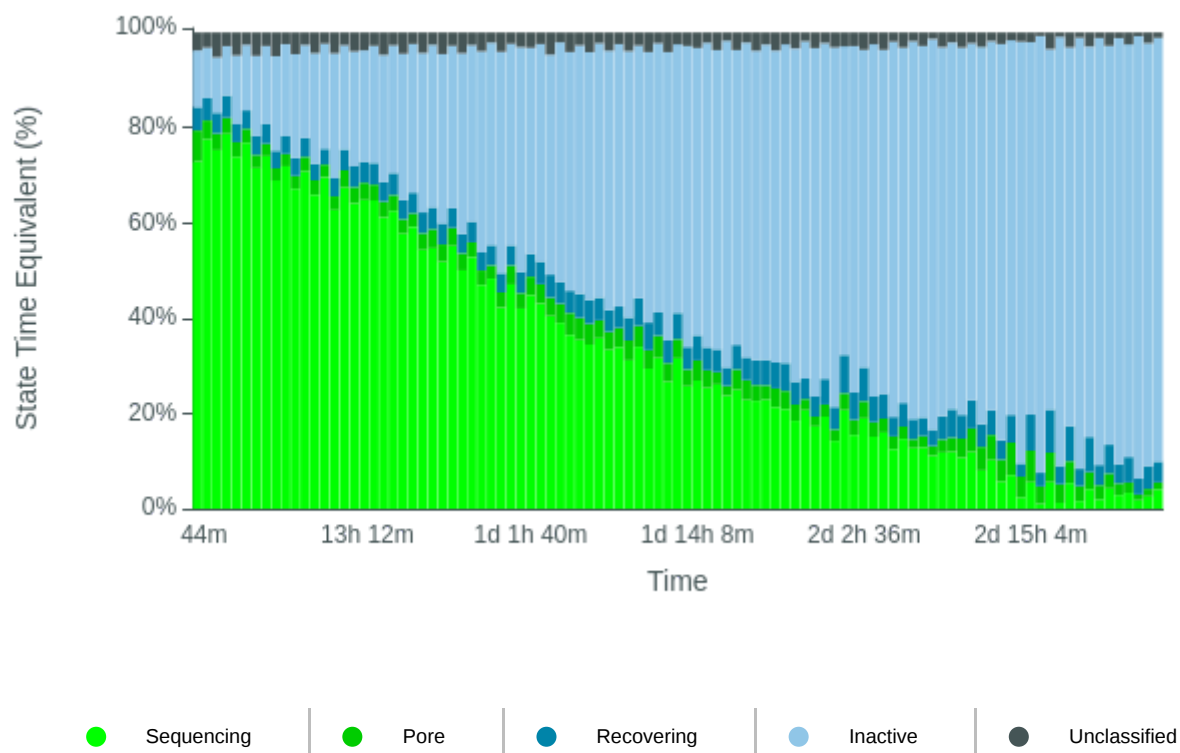


### Duty time Categorised

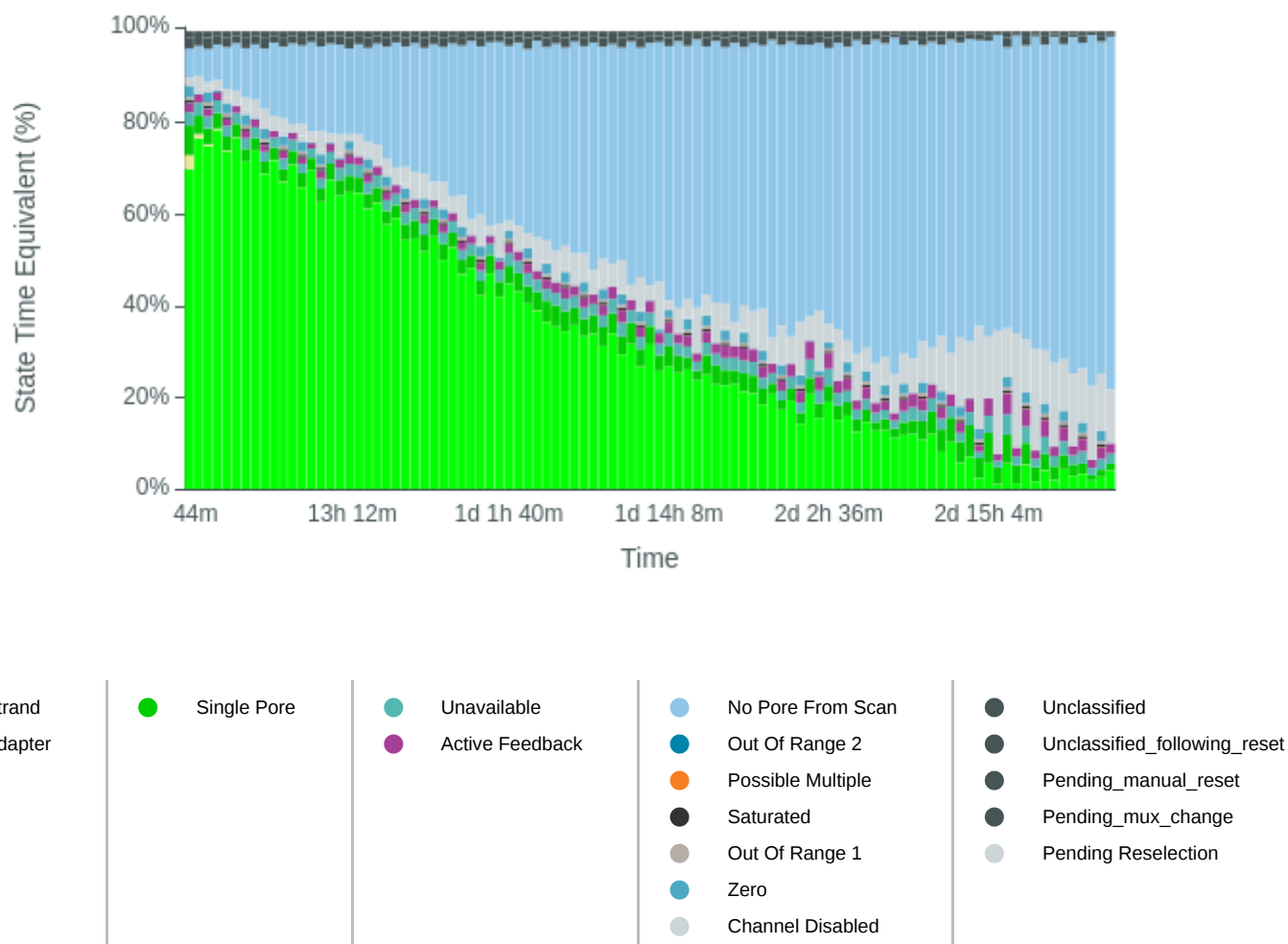


Mux Scan Grouped

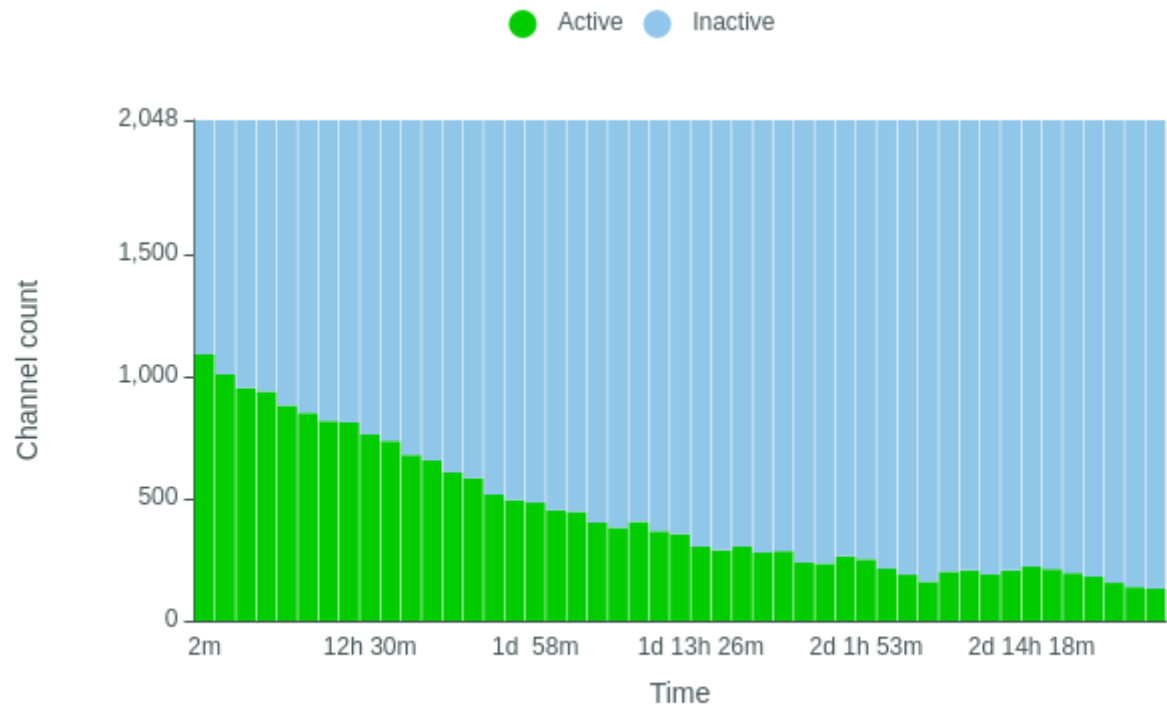


Mux Scan Categorised

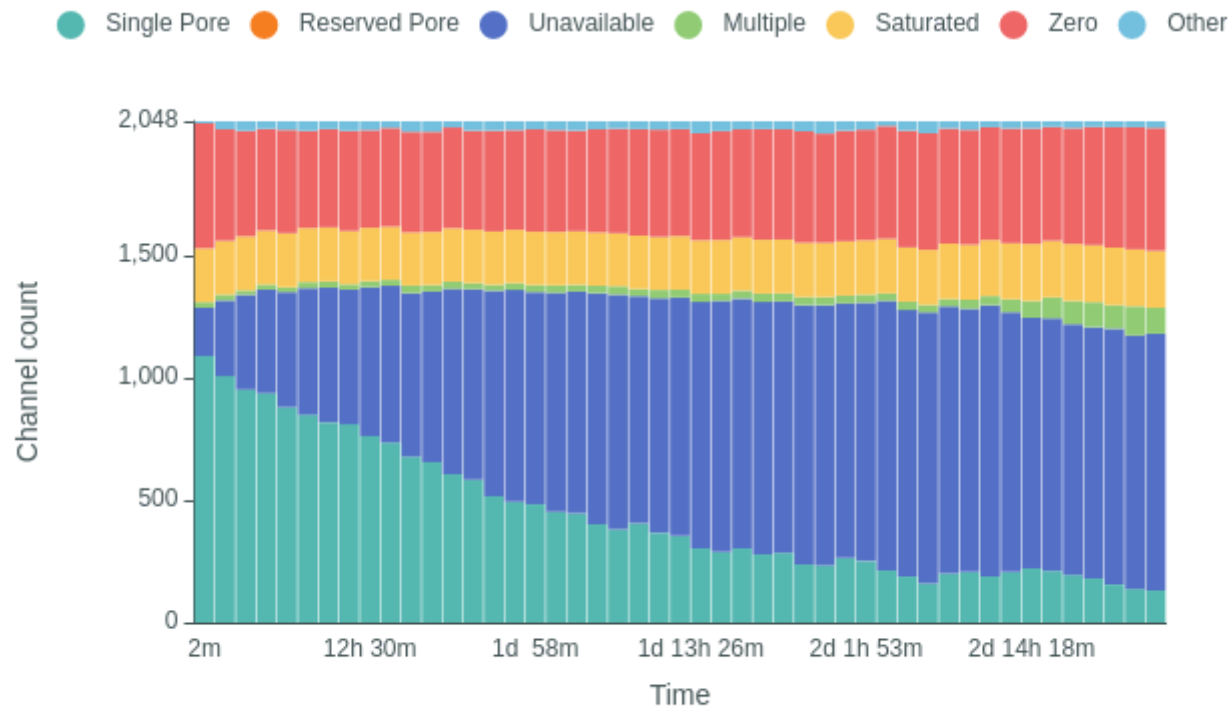


Temperature History

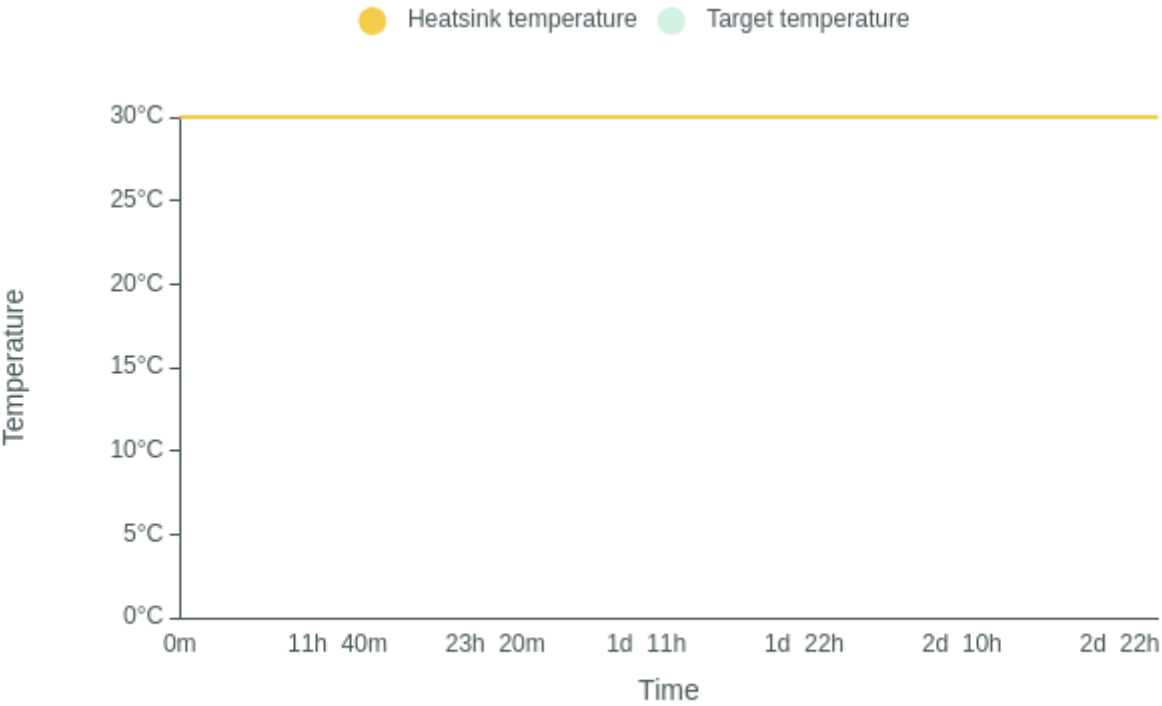


Bias Voltage History

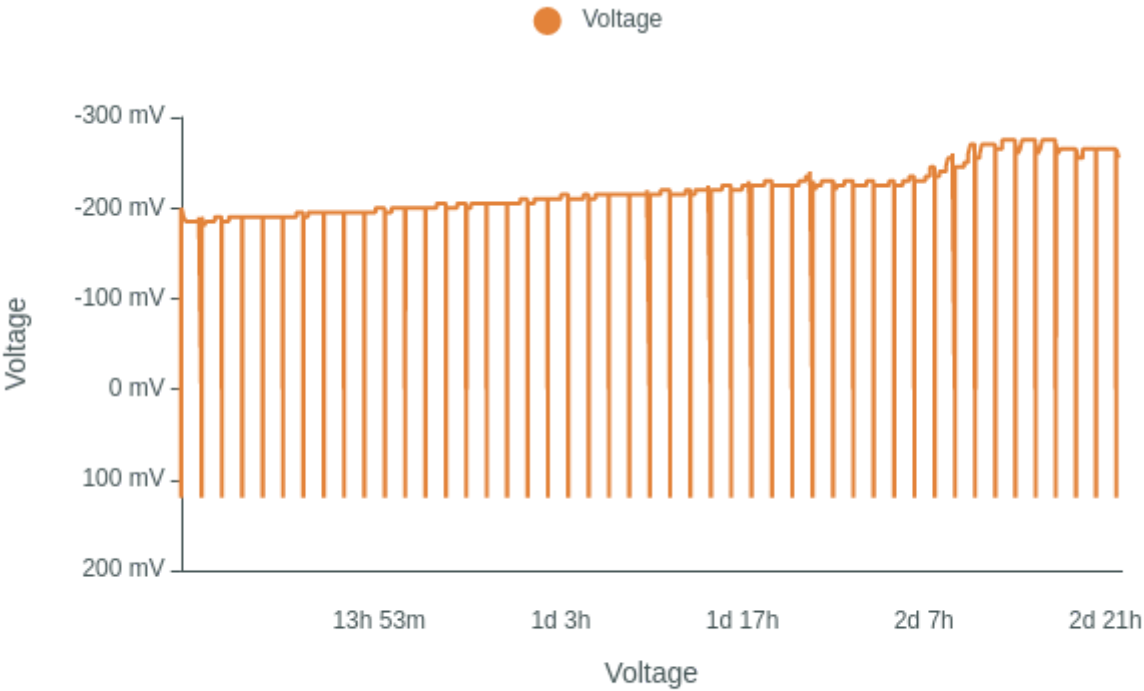

Translocation Speed

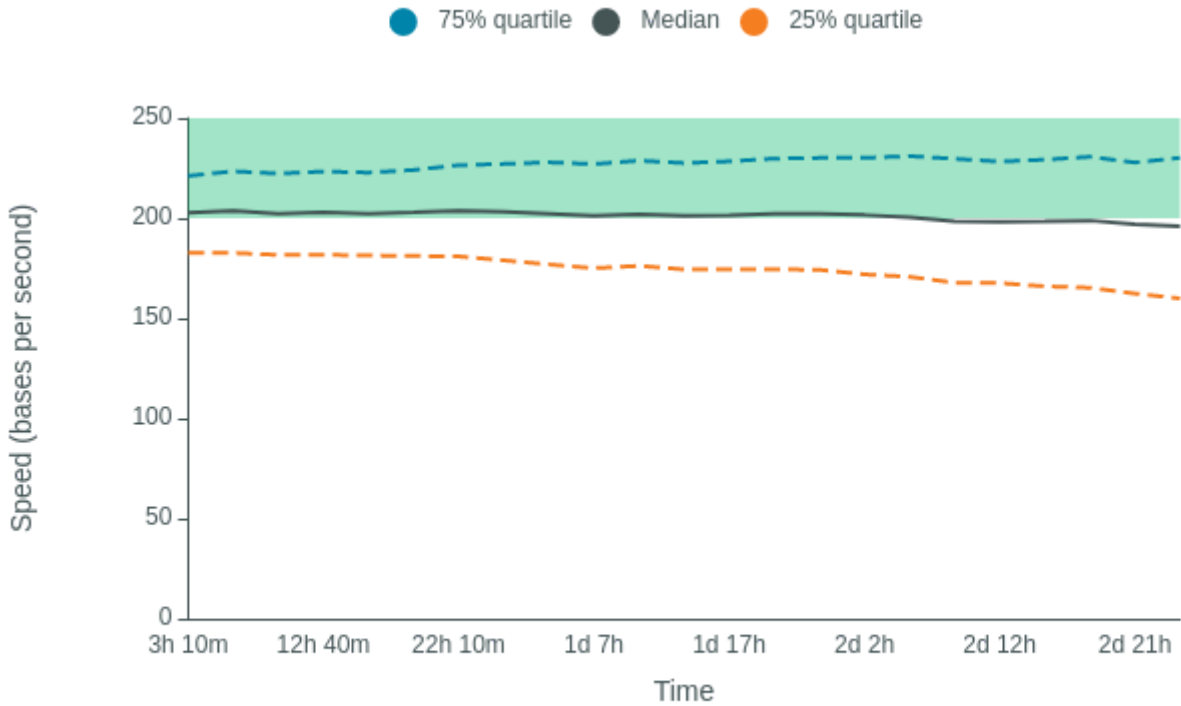


QScore

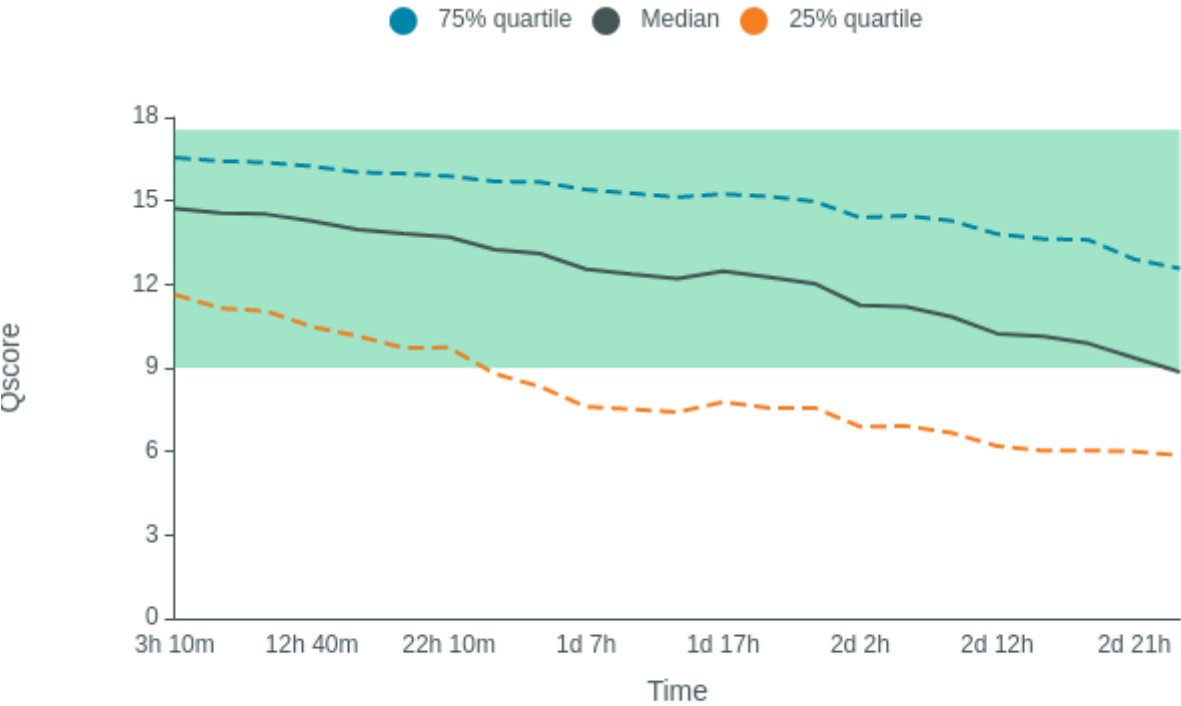



Disk Write Performance

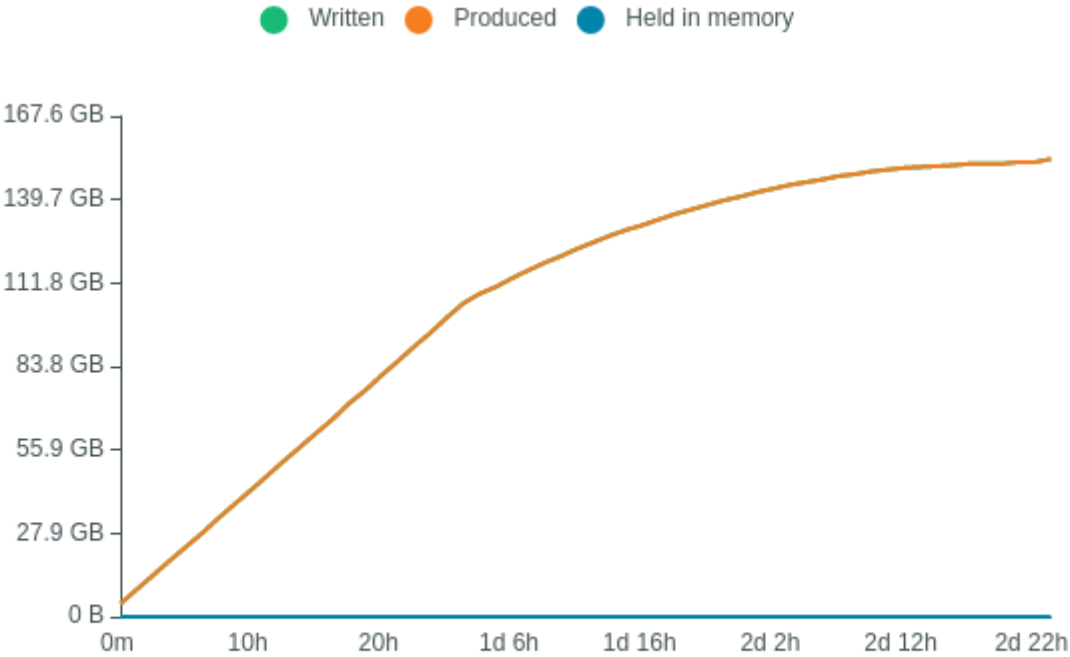

**Run Debug Messages**

- Mux scan for flow cell FAU08590 has found a total of 131 pores. 117 pores available for immediate sequencing June 13, 13:52
- Performing Mux Scan June 13, 13:50
- Mux scan for flow cell FAU08590 has found a total of 138 pores. 120 pores available for immediate sequencing June 13, 12:20
- Performing Mux Scan June 13, 12:17
- Mux scan for flow cell FAU08590 has found a total of 155 pores. 135 pores available for immediate sequencing June 13, 10:47
- Performing Mux Scan June 13, 10:45
- Mux scan for flow cell FAU08590 has found a total of 180 pores. 147 pores available for immediate sequencing June 13, 09:15
- Performing Mux Scan June 13, 09:12
- Mux scan for flow cell FAU08590 has found a total of 195 pores. 162 pores available for immediate sequencing June 13, 07:42
- Performing Mux Scan June 13, 07:40
- Mux scan for flow cell FAU08590 has found a total of 211 pores. 179 pores available for immediate sequencing June 13, 06:10
- Performing Mux Scan June 13, 06:07
- Mux scan for flow cell FAU08590 has found a total of 221 pores. 181 pores available for immediate sequencing June 13, 04:37
- Performing Mux Scan June 13, 04:35
- Mux scan for flow cell FAU08590 has found a total of 206 pores. 173 pores available for immediate sequencing June 13, 03:05
- Performing Mux Scan June 13, 03:02
- Mux scan for flow cell FAU08590 has found a total of 189 pores. 161 pores available for immediate sequencing June 13, 01:32
- Performing Mux Scan June 13, 01:30
- Mux scan for flow cell FAU08590 has found a total of 206 pores. 169 pores available for immediate sequencing June 12, 23:59
- Performing Mux Scan June 12, 23:57
- Mux scan for flow cell FAU08590 has found a total of 201 pores. 157 pores available for immediate sequencing June 12, 22:26
- Performing Mux Scan June 12, 22:24
- Mux scan for flow cell FAU08590 has found a total of 160 pores. 135 pores available for immediate sequencing June 12, 20:53
- Performing Mux Scan June 12, 20:50
- Mux scan for flow cell FAU08590 has found a total of 188 pores. 150 pores available for immediate sequencing June 12, 19:19
- Performing Mux Scan June 12, 19:17
- Mux scan for flow cell FAU08590 has found a total of 213 pores. 161 pores available for immediate sequencing June 12, 17:46
- Performing Mux Scan June 12, 17:43
- Mux scan for flow cell FAU08590 has found a total of 251 pores. 191 pores available for immediate sequencing June 12, 16:12
- Performing Mux Scan June 12, 16:10
- Mux scan for flow cell FAU08590 has found a total of 264 pores. 208 pores available for immediate sequencing June 12, 14:39
- Performing Mux Scan June 12, 14:37
- Mux scan for flow cell FAU08590 has found a total of 234 pores. 182 pores available for immediate sequencing June 12, 13:06

- Performing Mux Scan June 12, 13:03
- Mux scan for flow cell FAU08590 has found a total of 238 pores. 179 pores available for immediate sequencing June 12, 11:32
- Performing Mux Scan June 12, 11:30
- Mux scan for flow cell FAU08590 has found a total of 285 pores. 212 pores available for immediate sequencing June 12, 09:59
- Performing Mux Scan June 12, 09:56
- Mux scan for flow cell FAU08590 has found a total of 279 pores. 200 pores available for immediate sequencing June 12, 08:25
- Performing Mux Scan June 12, 08:23
- Mux scan for flow cell FAU08590 has found a total of 304 pores. 220 pores available for immediate sequencing June 12, 06:52
- Performing Mux Scan June 12, 06:49
- Mux scan for flow cell FAU08590 has found a total of 288 pores. 212 pores available for immediate sequencing June 12, 05:18
- Performing Mux Scan June 12, 05:16
- Mux scan for flow cell FAU08590 has found a total of 303 pores. 214 pores available for immediate sequencing June 12, 03:45
- Performing Mux Scan June 12, 03:42
- Mux scan for flow cell FAU08590 has found a total of 354 pores. 241 pores available for immediate sequencing June 12, 02:11
- Performing Mux Scan June 12, 02:09
- Mux scan for flow cell FAU08590 has found a total of 366 pores. 240 pores available for immediate sequencing June 12, 00:38
- Performing Mux Scan June 12, 00:35
- Mux scan for flow cell FAU08590 has found a total of 405 pores. 266 pores available for immediate sequencing June 11, 23:04
- Performing Mux Scan June 11, 23:02
- Mux scan for flow cell FAU08590 has found a total of 381 pores. 258 pores available for immediate sequencing June 11, 21:31
- Performing Mux Scan June 11, 21:28
- Mux scan for flow cell FAU08590 has found a total of 402 pores. 275 pores available for immediate sequencing June 11, 19:57
- Performing Mux Scan June 11, 19:55
- Mux scan for flow cell FAU08590 has found a total of 446 pores. 280 pores available for immediate sequencing June 11, 18:24
- Performing Mux Scan June 11, 18:21
- Mux scan for flow cell FAU08590 has found a total of 453 pores. 292 pores available for immediate sequencing June 11, 16:50
- Performing Mux Scan June 11, 16:48
- Mux scan for flow cell FAU08590 has found a total of 483 pores. 307 pores available for immediate sequencing June 11, 15:17
- Performing Mux Scan June 11, 15:14
- Mux scan for flow cell FAU08590 has found a total of 493 pores. 307 pores available for immediate sequencing June 11, 13:43
- Performing Mux Scan June 11, 13:41
- Mux scan for flow cell FAU08590 has found a total of 517 pores. 313 pores available for immediate sequencing June 11, 12:10
- Performing Mux Scan June 11, 12:07
- Mux scan for flow cell FAU08590 has found a total of 584 pores. 341 pores available for immediate sequencing June 11, 10:36
- Performing Mux Scan June 11, 10:34

- Mux scan for flow cell FAU08590 has found a total of 607 pores. 358 pores available for immediate sequencing June 11, 09:03
- Performing Mux Scan June 11, 09:01
- Mux scan for flow cell FAU08590 has found a total of 656 pores. 368 pores available for immediate sequencing June 11, 07:29
- Performing Mux Scan June 11, 07:27
- Mux scan for flow cell FAU08590 has found a total of 678 pores. 372 pores available for immediate sequencing June 11, 05:56
- Performing Mux Scan June 11, 05:54
- Mux scan for flow cell FAU08590 has found a total of 735 pores. 399 pores available for immediate sequencing June 11, 04:22
- Performing Mux Scan June 11, 04:20
- Mux scan for flow cell FAU08590 has found a total of 762 pores. 411 pores available for immediate sequencing June 11, 02:49
- Performing Mux Scan June 11, 02:47
- Mux scan for flow cell FAU08590 has found a total of 810 pores. 411 pores available for immediate sequencing June 11, 01:16
- Performing Mux Scan June 11, 01:13
- Mux scan for flow cell FAU08590 has found a total of 818 pores. 412 pores available for immediate sequencing June 10, 23:42
- Performing Mux Scan June 10, 23:40
- Mux scan for flow cell FAU08590 has found a total of 849 pores. 421 pores available for immediate sequencing June 10, 22:09
- Performing Mux Scan June 10, 22:06
- Mux scan for flow cell FAU08590 has found a total of 880 pores. 430 pores available for immediate sequencing June 10, 20:35
- Performing Mux Scan June 10, 20:33
- Mux scan for flow cell FAU08590 has found a total of 937 pores. 450 pores available for immediate sequencing June 10, 19:02
- Performing Mux Scan June 10, 18:59
- Mux scan for flow cell FAU08590 has found a total of 952 pores. 459 pores available for immediate sequencing June 10, 17:28
- Performing Mux Scan June 10, 17:26
- Mux scan for flow cell FAU08590 has found a total of 1007 pores. 472 pores available for immediate sequencing June 10, 15:55
- Performing Mux Scan June 10, 15:52
- Mux scan for flow cell FAU08590 has found a total of 1089 pores. 478 pores available for immediate sequencing June 10, 14:21
- Performing Mux Scan June 10, 14:19
- Starting sequencing procedure June 10, 14:19
- Waiting up to 300 seconds for temperature to stabilise at 30.0°C June 10, 14:15